BIOS living

ROLLATOR
56045

Instruction
Manual

BIOS living

The BIOS Living rollator comes with removable nylon storage bags, locking brakes, and a padded seat with back rest. It folds away easily for storage and transport and comes with anti-fold safety locks to prevent the rollator from accidentally folding, while in use.

Please read these instructions before assembly or use of this rollator, and please keep the manual in a safe place.

DO NOT assemble or use this equipment without first reading and understanding this instruction manual.

Personal Safety

- The rollator can provide assistance to an individual weighing up to 136 kg (300 lbs), including the weight of the contents of the nylon storage bags
- The maximum load of the nylon basket under seat is 10 kg (22 lbs)
- Never stand on the rollator, or use it as a step
- Never transport a person on the seat of the rollator
- Use caution when lifting and carrying the rollator
- Test the brakes before using. Care should be taken to ensure that ALL hand and height adjustments are secured, and that moving parts are in good working order before using this or any mobility aid
- Ensure that brakes are in the locked position before sitting
- Never sit on the rollator when it is on an incline
- Always sit in the center of the seat
- Never sit on the rollator when it is folded
- Each individual should always consult with their physician or

therapist to determine proper adjustment usage


- Avoid excessive weight in the rollator bag, it may cause rollator to tip
- Max weight load for rollator bag: 0.91 kg (2 lbs)
- Customers should seek the advice / recommendations of a medical professional, to determine if the rollator is suitable for their needs.

Assembly Instructions



Please examine the rollator for shipping damage. Do not use the rollator if there is apparent damage. Please contact your place of purchase for repair and/or return information.

1. Before assembling read all instructions thoroughly and make sure all the parts listed below are included.


- A. Right Handle
- B. Left Handle
- C. Backrest
- D. Front Wheels (2)
- E. Rollator Frame
- F. Rear Wheels (2)
- G. Fastening Knobs (2)
- H. Safety Buckles (2)
- I. Nylon Basket
- J. Rollator Bag
- K. Stability Knob (2)

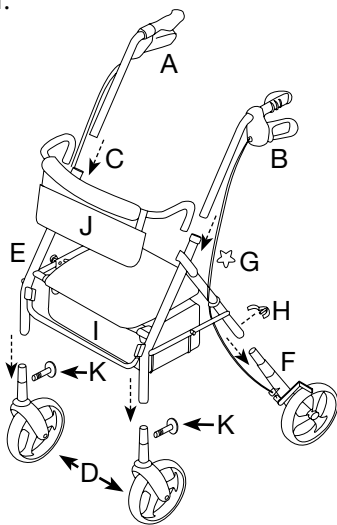
2. Expand the rollator frame by pulling outwards on the frame lock. Push the BIOS Safety Lock down and into place to lock the cross bars. To fold the rollator pull up on the BIOS Safety Lock to disengage. Place the rollator frame on a flat, hard surface. Remove plastic wrap on the wheels.

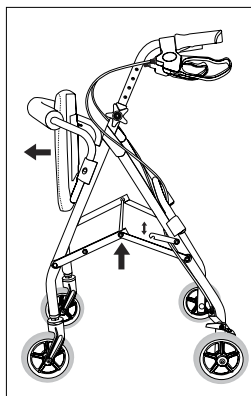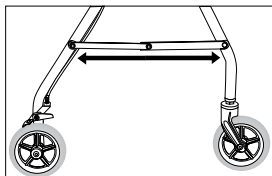
3. Insert handles (A and B) into the frame as illustrated. Insert the triangular fastening screw through the frame and handle holes, use the triangular fastening knobs (G) to adjust the handle height that best suits you; ideally the handles should be positioned at wrist height. Secure the handles in place by turning the knobs clockwise until tightened (do not over tighten).

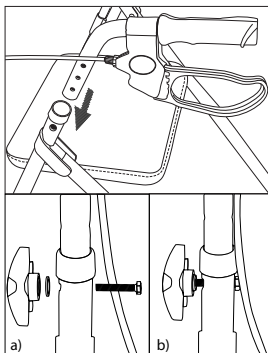4. Attach the rear wheels (F) into the frame, ensure that the metal beads are aligned correctly with the holes. Clip safety buckles (H) onto rear legs. Attach the front wheels (D) in the same manner. (K) Twist stability knobs into holes located at the rear of the front legs.

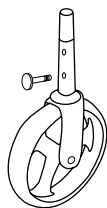
5. Insert the curved backrest (C) into its slot on the front side of the rollator frame as illustrated.

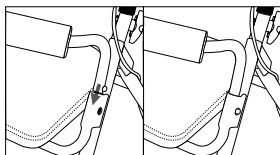6. Attach the nylon rollator bag to the backrest (optional).


Brake Operation

To test the brakes, hold both brakes and gently squeeze the brake handle. The brakes should activate when this is done. If adjustments are required, turn the screw and nuts located at the end of the brake cable to increase (counter clockwise) or decrease (clockwise) the brake tension.

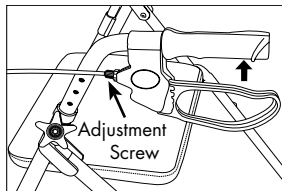
Folding the Rollator

Lift the seat, and pull up on the metal bar which is attached to the cross folding bar and pull upwards on the safety lock to release. The nylon storage bag will fold up into the rollator.

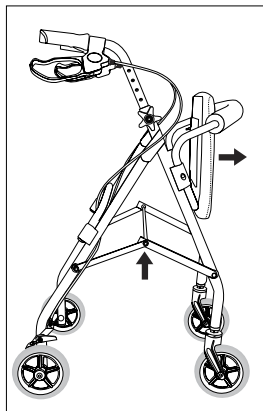Parking Brakes

1. Push both handles down until they click to activate parking brake. The handle will lock automatically and stay in place.
2. To release the parking brake, pull the handle up. The handle will return to its basic position.

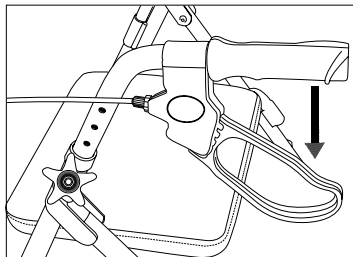Parking Brake

Care & Maintenance

1. Ensure that ALL attached hardware is secured at all times.
2. Make sure seat is secured to rollator properly, tighten screw if loose.
3. Replace any broken, damaged or worn items immediately.

Recommended Handle/Seat Height Adjustment

Handle Height	Approx. User Height Min/Max
79.5 - 95 cm (31.3" - 37.4")	162 - 176 cm (5'-4" - 5'-10")
Seat Height	
49 - 53 cm (19.29" - 20.86")	

One Year Warranty

If this product proves to be defective in material or workmanship within one year of purchase, please return it to the address below. It will be repaired or replaced without charge upon receipt of the unit prepaid with \$20.00 to cover handling, packaging and return postage. Please include proof of purchase, your full name, address, daytime phone number or email address.

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